


Polyamide (15) Nylon 6/10 (HK-404L)

Nylon 6/10 was analyzed using GPC HK-404L, an organic SEC (GPC) column for ultra-rapid analysis.

Sample : Nylon 6/10 (Polyhexamethylene sebacamide) 0.2 %, 5 μ L

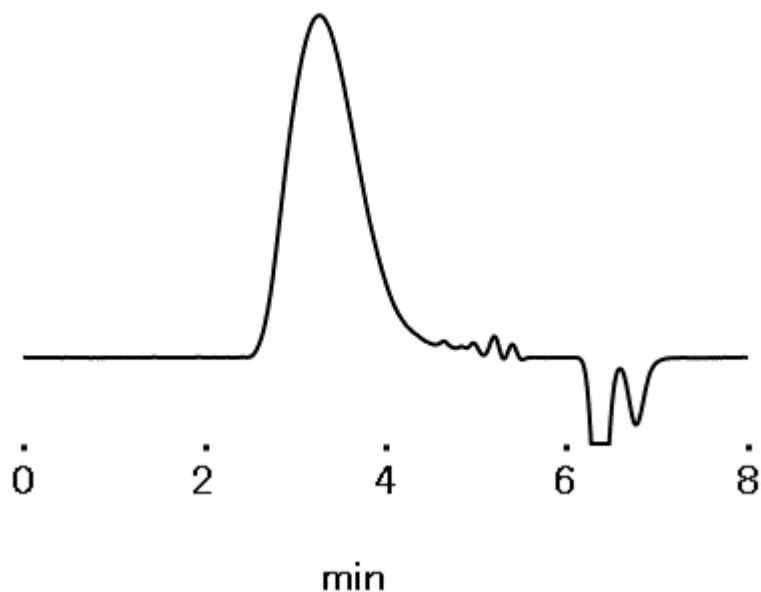

Column	: Shodex GPC HK-404L (4.6 mm I.D. x 150 mm)
Eluent	: 5 mM CF ₃ COONa in HFIP
Flow rate	: 0.3 mL/min
Detector	: RI (small cell volume)
Column temp.	: 40 °C