


Polyamide (16) Nylon 6/6 (HK-404L)

Nylon 6/6 was analyzed using GPC HK-404L, an organic SEC (GPC) column for ultra-rapid analysis.

Sample : 5 μ L

Nylon6/6

(Polyhexamethylene adipamide) 0.22 %

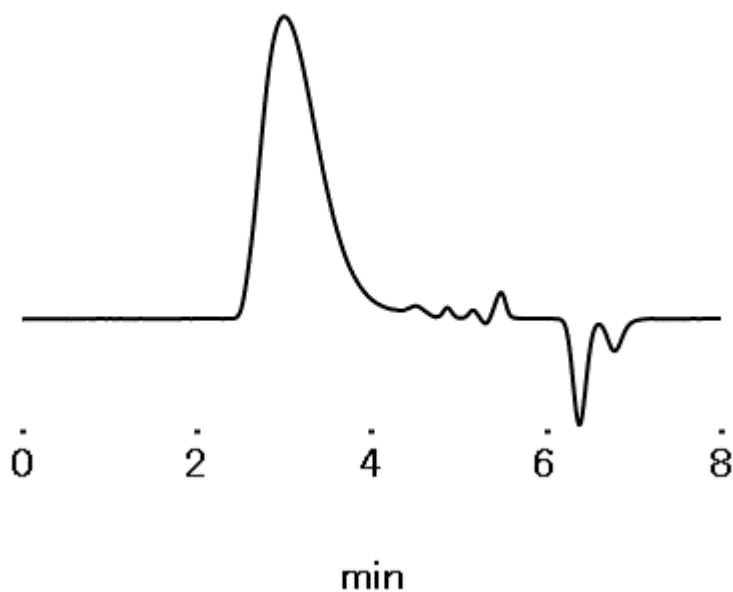

Column	: Shodex GPC HK-404L (4.6 mm I.D. x 150 mm)
Eluent	: 5 mM CF ₃ COONa in HFIP
Flow rate	: 0.3 mL/min
Detector	: RI (small cell volume)
Column temp.	: 40 °C