


Polyamide (18) Nylon 6(3)T (HK-404L)

Nylon 6(3)T was analyzed using GPC HK-404L, an organic SEC (GPC) column for ultra-rapid analysis.

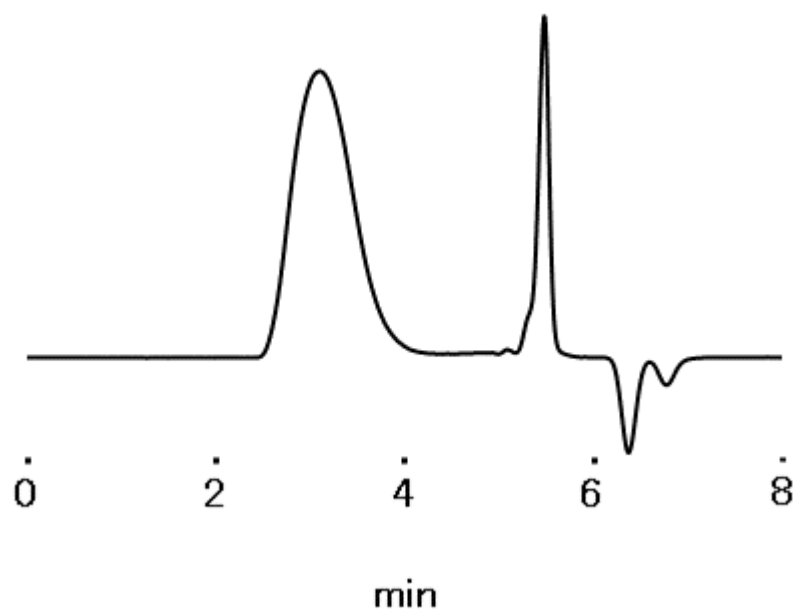

Sample : 5 μ L
Nylon6(3)T 0.2 %

Column	: Shodex GPC HK-404L (4.6 mm I.D. x 150 mm)
Eluent	: 5 mM CF_3COONa in HFIP
Flow rate	: 0.3 mL/min
Detector	: RI (small cell volume)
Column temp.	: 40 °C